HPLC Application

Distearin in Palm Oil on Phenogel 5u 100A


Column:	Phenogel™ 5 µm 100A, LC Column 300 x 7.8 mm, Ea		
Dimensions:	300 x 7.8 mm ID		
Order No:	00H-0442-K0		
Elution Type:	Isocratic		
Eluent A:	THF		
Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100
Flow Rate:	1 mL/min		
Col. Temp.:	ambient		
Detection: 8330	Differential refractometer @ nm (ambient)		



Ophenomenex[®]

Products used in this application:



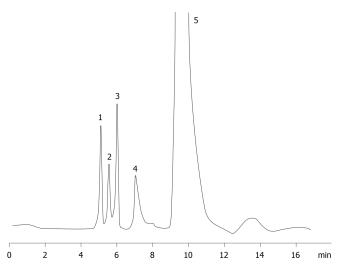

ANALYTES:

- 1 Palm oil
- 2 Distearin
- Palm oil methyl ester 3
- 4 Monostearin
- Methanol 5

©2025 Phenomenex Inc. All rights reserved.

Phenomenex products are available worldwide.

www.phenomenex.com.cn

For more information contact your Phenomenex Representative at info@phenomenex.com